




MASSAGE CHAIR USER MANUAL

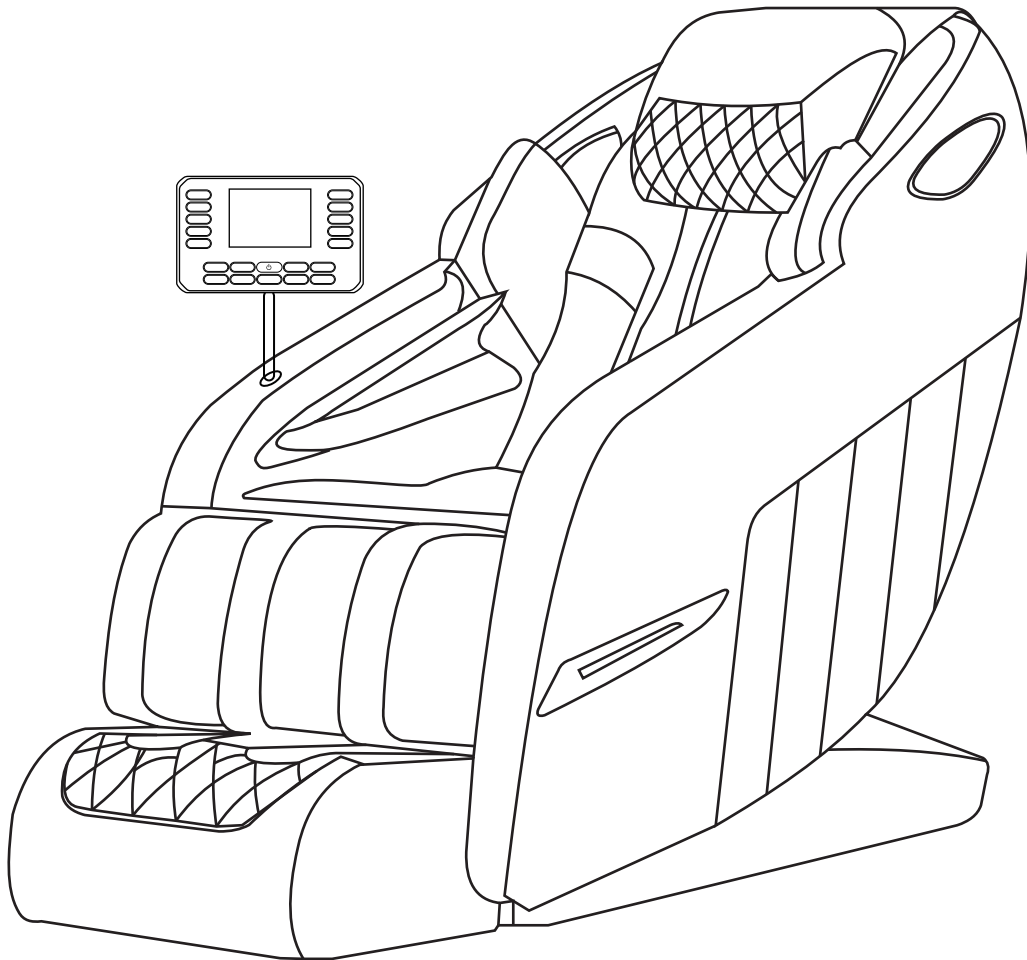

Please read SAFETY INSTRUCTIONS carefully before using the chair.

CONTENTS

Please read and follow all safety and operating instructions listed in this user manual prior to use.
Basic precautions should always be used when operating.

Safety Instructions	- - - - -	3-4
Product Specifications	- - - - -	5
Product Components	- - - - -	6
Before Use	- - - - -	7-8
Instructions	- - - - -	9-11

Note: Our company reserves the right to revise the design and description of this product without any further notice.

SAFETY INSTRUCTIONS

Warning

The plug must be fully inserted into an appropriate outlet. Improper connection of the equipment-grounding conductor can result in the risk of electric shock

Check with a qualified electrician or service technician, if you are in doubt as to whether The product is properly grounded. Do not modify the plug provided with the product. If the outlet is not suitable. Please ensure that a proper outlet is installed by a qualified electrician.

When using the massage chair, basis precautions should always be followed to reduce the risk of electric shock, burns, fire, or injury to person.

1. Please ensure that you plug the massage chair into a proper electrical outlet.
2. Make sure the power plug is securely inserted into the proper electrical outlet to reduce the risk of short-circuit and fire.
3. Always unplug the massage chair from the electrical outlet immediately after use and before cleaning.
4. Do not use the massage chair in places where the humidity is high.
5. Do not use the massage chair outdoors and avoid place with direct sunlight.
6. If the power goes off during use, be sure to turn the power switch to the off position and pull the power plug out from the socket to prevent damage from occurring when the power is restored.
7. Never operate the massage chair if the power cord is wet or damaged.
8. Never plug/unplug the power cord from the socket with wet hands.
9. Close and continuous supervision is necessary when this appliance is being used by on or near children invalid or disabled person.
10. Use this appliance only for its intended use as described in this use manual.
11. Do not use any parts not from factory.
12. Do not stand, sit or jump in the calves rest when massage chair is in use to prevent any injury to person.
13. Do not attempt to stand on the backrest and armrest to prevent any injury to person.
14. Be carefully not to fall asleep during the massage session.
15. Always keep the massage chair clean and never insert any object into the gaps on and Around the massage chair.
16. Do not move the massage chair by pulling the power supply cord.
17. Always check whether the surface of the massage chair is torn or ripped before usage, if you find any unfamiliar rips in any area of the cloth in the backrest, do not use the massage chair, kindly contact the service centre immediately for required.
18. This appliance is designed for household use only.
19. For those suffering from illnesses, please consult your doctor for advice prior to use.

20.The following persons are requested to consult a doctor prior to use in order to reduce the risk of user feeling sick or suffering from injury

- Those who have malignant tumor
- . Those who suffer from heart disease
- .Those who are experiencing pregnancy, menstruation and puerperal fever
- . Those who require complete rest
- .Those who are under medical treatment
- Those who suffer from back problem or whose backbone is abnormal or bent

21. If you experience pain during usage, stop using the massage chair immediately, if any injuries occur, seek medical treatment immediately.

22. Never massage the swollen or inflamed parts of the body.

23. This product is non-professional appliance designed to provide a comfortable massage. do not use this product as substitute for medical treatment.

24. Each massage session should not exceed 30 minutes. massage the targeted area only for 15mins, over-massaging can lead to overheating , thus giving the appliance a shorter lifespan.

25. Do not operate where aerosol(spray) product are used or where oxygen is being administered.

26.The massage chair is not suitable for children aged 12 and below.

NOTE:

·This appliance is not intended for use by persons(including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

.If the supply cord is damaged,it must be replaced by manufacturer,it's service agent or similarly qualified persons in order to avoid a hazard.


·The appliance heats up during operation. persons sensitive to heat must take extra caution when using the appliance.

Product Specifications


Description	Specification
Product	Massage chair
Model	CALMNEST
Voltage	110-240V
Power	150W
Rated frequency	50-60HZ
Rated working time	15Mins

PRODUCT COMPONENTS

Main Body Diagram(Front)

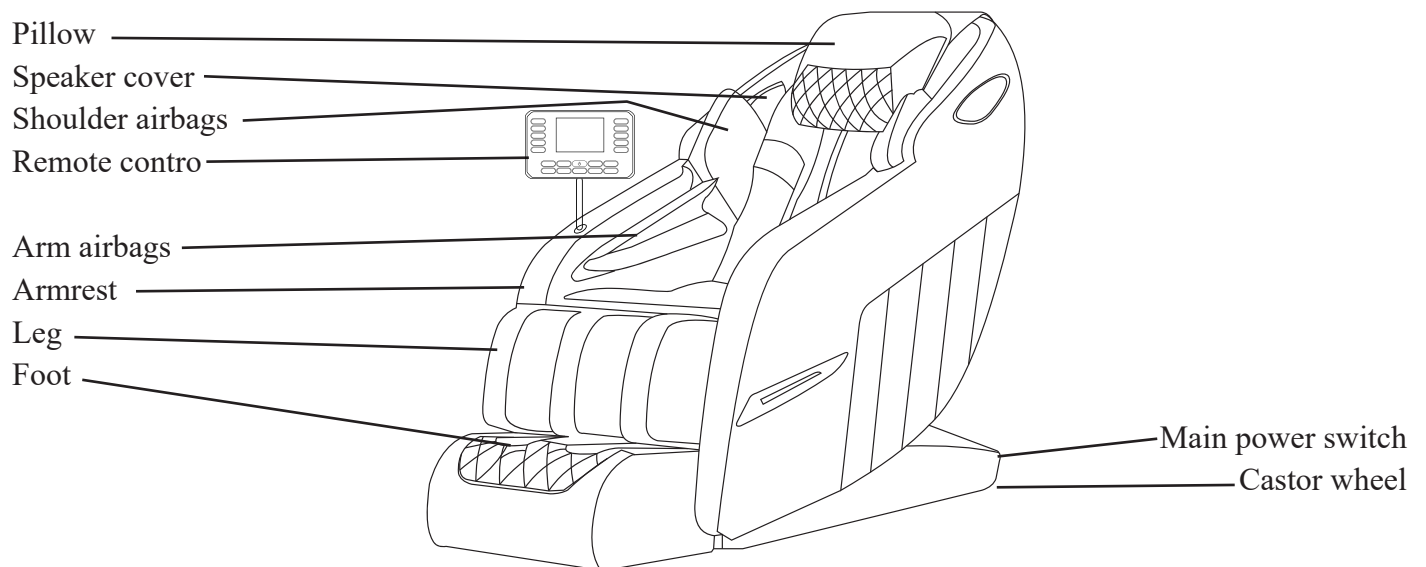


Main Body Diagram(Back)

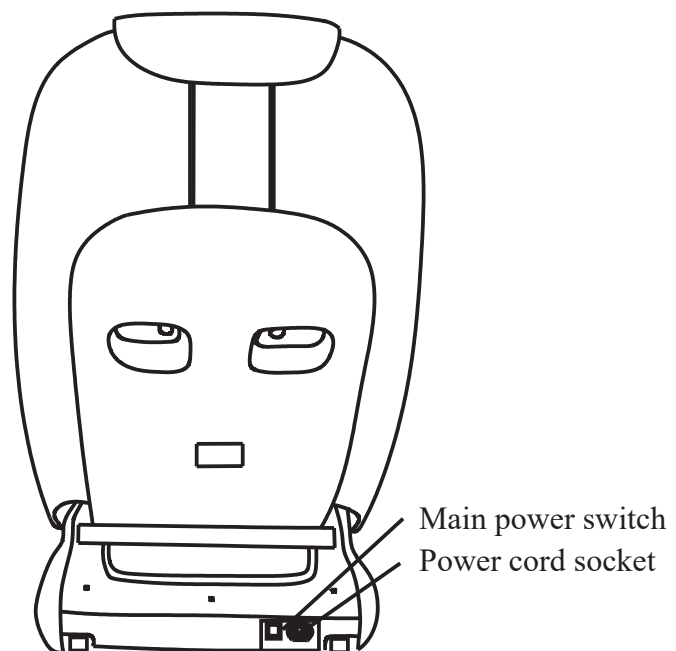

Before use

1.Remove the package and get out the spare parts

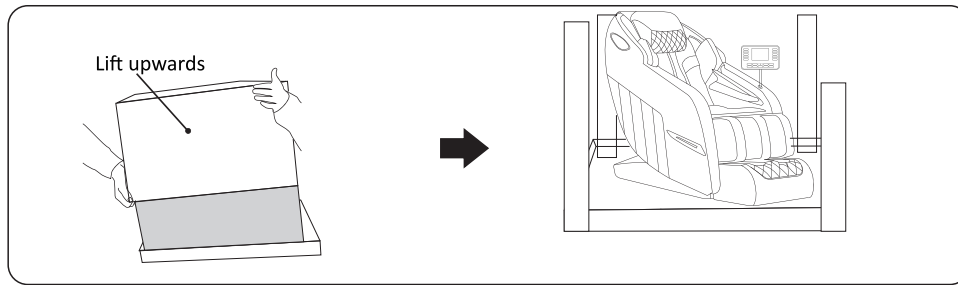

1)Cut off the packing strap and lift up and move away the carton.

2) Take out the spare parts box for next use.

Note:Please gently lower your legs from your hips to the ground.

2.Place the massage chair in an open area


There should be at least 25cm before chair and at least 50cm behind chair when place the chair, so the chair can be recline freely.

*In order to avoid the damage of the floor and to reduce the noise, please put area rug under the massage chair.

*Make sure that there are no barriers in front and the back of the massage chair.

*Please do not put it in wet or uneven floor.

*To avoid signal interference, please place it 4 feet away from AV,TV, and radio.

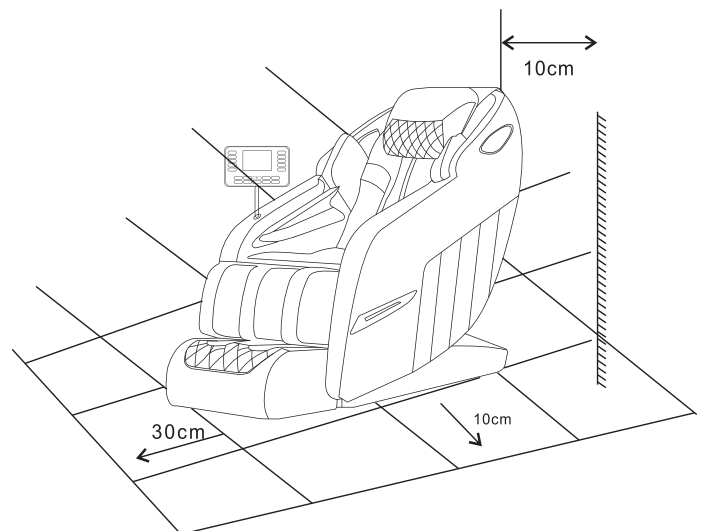

Warning:

Do not put the power cord under of relax chair!

3. Moving the massage chair.


As illustrated, push the backrest downwards and tilt chair on to its wheels at a 45°angle.

Slowly push and maneuver the massage chair to desired location. Lower massage chair slowly and smoothly until it is fully flat on the ground.

CAUTION:

*Do not move massage chair when occupied or in use.

*Moving the massage chair may cause damage to the floors. Floor protection is recommended (e.g. rugs or carpet).

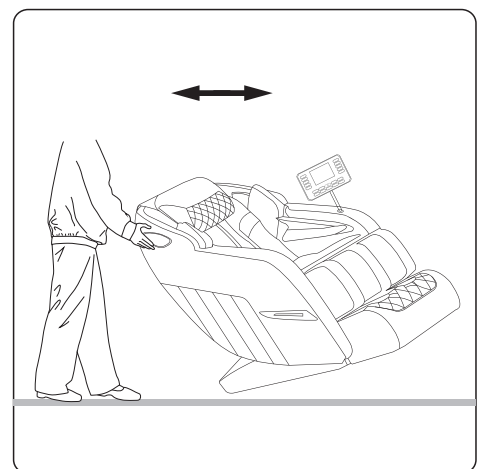


*Do not drop or suddenly let go of the chair when moving as it may damage the chair.

*Two persons are recommended when moving and extra care when handling the massage chair.

Before use


1. Make sure the product is in good condition and proper working order.

- Check the headrest and backrest pad.
- Inspect the power plug.
- Make sure the footrest is in its upright original position in order to prevent damage.

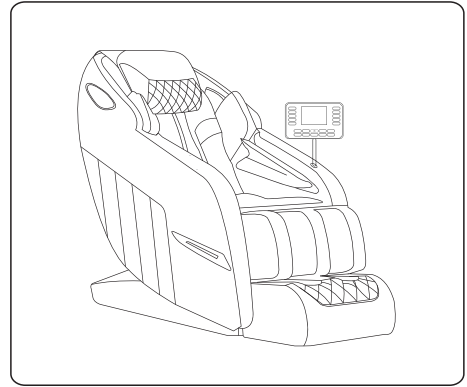


2. Check the surrounding environment

- Ensure there are no people, pets, or objects obstructing the massage chair.

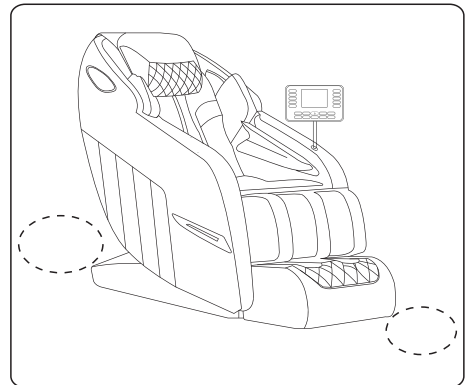


3. Plug in the massage chair power cord.

- Connect the power cord.
- Turn ON the power switch, located on the right side of the massage chair.

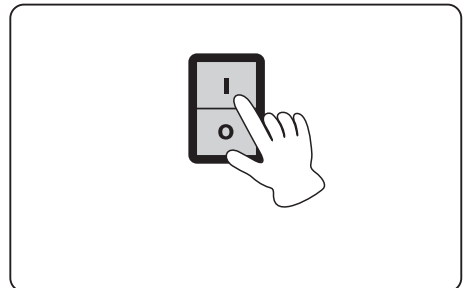


4. Adjust the position of the headrest.

- Adjust the headrest until it is in a comfortable position.

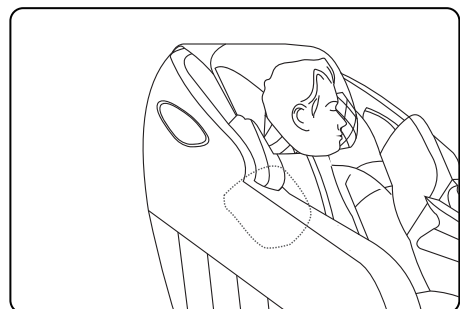

Instruction

Remote control

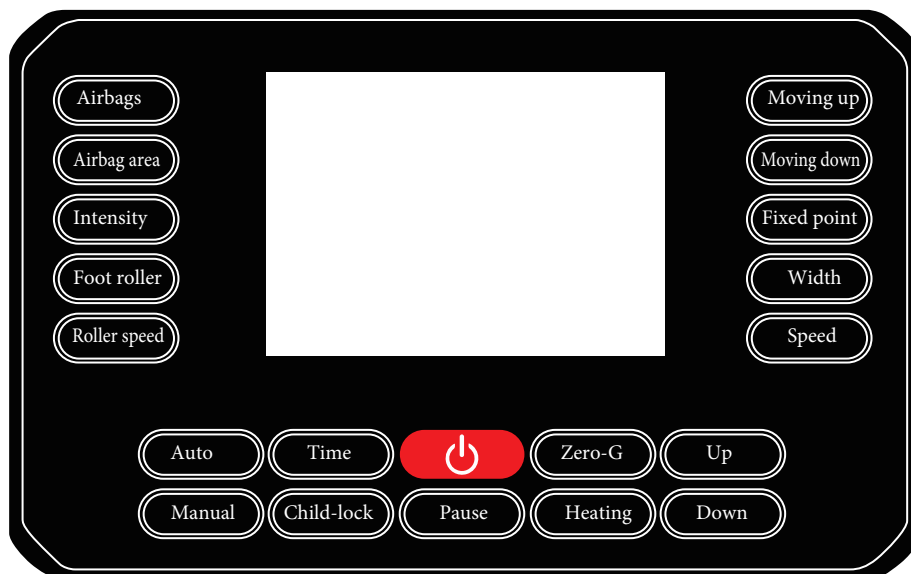

Press the "power" button to enter the automatic massage mode, and then press "Power" to stop the massage chair.

Airbags:

Press the "Airbags" button when the device is powered on to turn on or off the full-body airbags.

Airbag area:

Press the "Airbag area" button when the device is powered on to open airbags in all parts of the body.

Intensity:

Press the button to adjust the airbag intensity (1-3).

Foot roller:

Press the "Foot roller" button when the device is powered on to turn on the roller massage function; press it again to turn off the roller massage function.

Roller speed:

Press the "Roller speed" button when the machine is powered on to adjust the roller speed.

Auto:

By pressing the "Auto" button: you can choose a variety of different automatic massage programs to meet your different needs.

1. Relaxation
2. Shoulder peri-arthritis improvement
3. Back care
4. Stretching
5. Sports recovery
6. Powerful massage
7. Swing mode
8. Physical recovery
9. Sleep mode
10. Elderly mode
11. President mode
12. Queen mode
13. Thai massage
14. Chinese health care
15. Ancient massage
16. Full airbag.

Manual:

Press the "Manual" key when the machine is turned on. Turn on specific massage functions such as kneading, beating, kneading, patting, and acupressure.

Time:

Press the "Time" button to select different massage times when the machine is powered on.

Child-lock:

When the device is powered on, long press the "Child Lock" button to lock all buttons on the screen to prevent them from being touched. Long press again to unlock.

Pause:

Press the "Pause" button when the machine is on; the machine will stop working.

Zero-G:

Press the "Zero-G" button when the device is powered on to start the zero-gravity massage mode.

Heating:

Press this button to turn on or off the heating massage.

Up:

Press the "up" button under power-on condition, the back of the chair will be lifted up.

Down:

Press the "Down" button in the power-on state and recline the back of the chair.

Speed:

Press the "Speed" button when the machine is powered on to adjust the movement speed.

Width:

Press the "Width" button in manual mode to adjust the width, medium and narrowness of the distance between the massage heads.

Fixed point:

In manual mode, press the "Fixed Point" button to start fixed point massage.

Moving down:

Press the "Moving down" button in manual mode when the machine is powered on to adjust the movement to move downward.

Moving up:

Press the "Moving up" button in manual mode when the machine is powered on to adjust the movement to move upward.

Instruction

Pairing massage chair with bluetooth music

1).Pairing: turn on bluetooth. device in the handphone or ipad.

Search for massage chair.

2).Playing music: after pairing use the control device in the handphone or ipad to manage.

Note:disconnect Bluetooth. when the massage chair is not in use.

After massage

1.After massage press the on/off button to return the chair to its original position.



*Ensure nothing is obstructing the chair from returning to its original position.

2. Press the power switch of the control panel to shut down all the functions.

3. Press the power switch to off to cut off power supply.

4.Unplug from the wall socket.

* Power switch located the rear of the massage chair

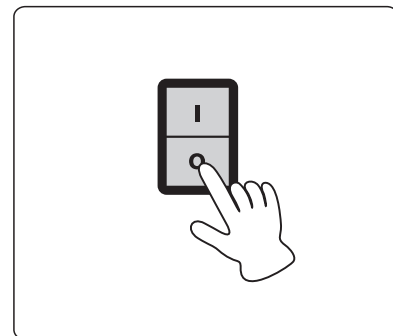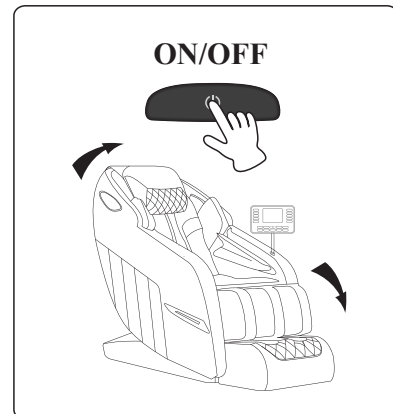

Trouble Shooting

In case of malfunction or emergence of abnormalities, stop using the massage chair and unplug from the power source. Try the following methods to resolve the problem. If the problem persists contact the service agent immediately.

Fault phenomenon	Possible reasons	Solution
There is a sound (operating sound) *The sound of the airbag inflating *The sound of the airbag deflation *The sound of the solenoid valve switching *The sound of the air pump working	The product makes mechanical sounds when working	This sound is normal and no action is required.
The whole machine does not work at all	The power plug is not plugged into the socket The socket switch is not turned on The massage button is not pressed	Check the power plug and socket, and make sure the plug and socket are in good contact. Press the massage button
The massage is in progress suddenly stopped	This product has just completed the rated working time of 15 minutes and automatically switched to standby mode	If you want to massage again, press the power button again to turn on the massage function
	The rated use time has been exceeded, and the motor is overheated and protected	Turn off the machine, cut off the power, and wait for the motor to cool before using it

MAINTENANCE

- Turn off and pull out the plug before clean the massage chair, do not touch the plug with wet hands.
- Use dry cloth to wipe sweat and other besmirch.
- Besmirch is serious, please use neutral wash to wipe and lightly brush the massage chair.
- Please airing.
- Prohibit to use water to clean the chair.
- Do not dry the chair with roaster.
- Prohibit to use alcohol, benzene wash etc.